


The most reliable VOC concentrator on the market today



Key Features


- Superior moisture control (patented feature)
- Low carryover (patented feature)
- Easy maintenance & diagnostics
- Best in-class service & support
- 3 year warranty on electronic boards

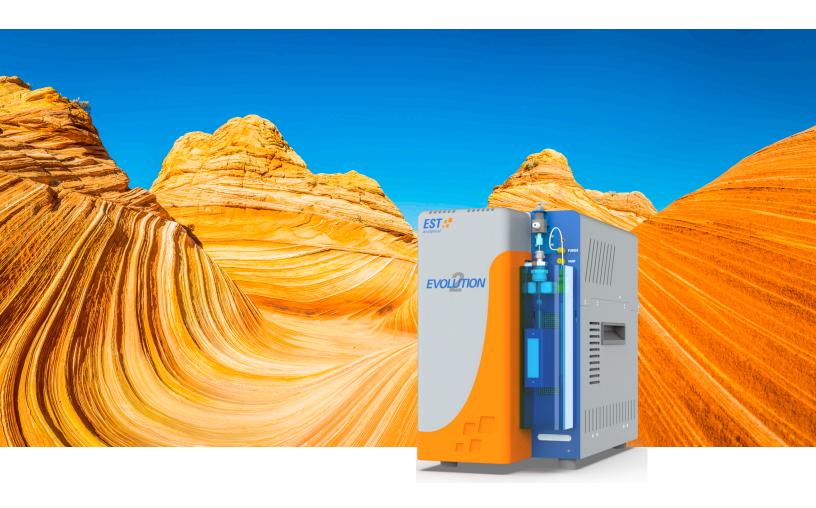

Methods

502.1 502.2	5030 5035
002.2	0000
524.2	8010
524.3	8015
524.4	8020
	8021
601	8030
602	0000
603	8240
	8260
624	0_00

Water Management

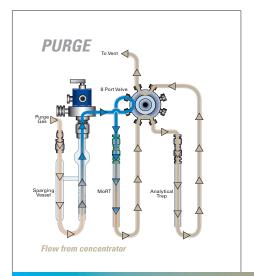




- The Evolution is designed to remove water during the purge process.
- By using an 8-port valve instead of a 6-port valve, the sample during desorption goes directly to the GC.



Water Management System

Superior Chromatography

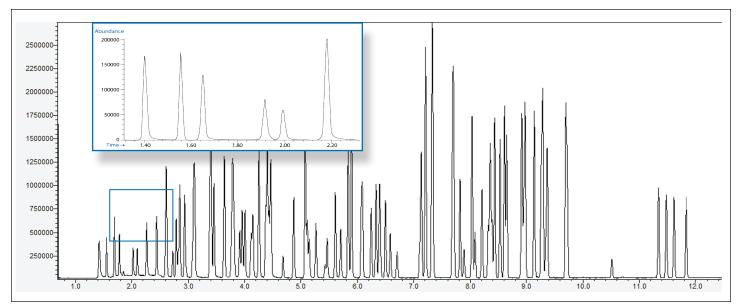

50 ppb EPA Method 8260. Front end gases showing improved resolution are highlighted.

Reduced Carryover



The eliminator is EST analytical's patented sparge heater that bakes the glassware of the Evolution concentrator between sample runs.







Visit estanalytical.com